

IncuWork *mi*



ADDRESS
28/9 MOO 1 BANGSRIMUANG,
MUANG, NONTHABURI 11000
THAILAND

TELEPHONE & FAX
+66 2886 7808
+66 2886 7809


E-MAIL
SALES@PULSESCIENCE.CO.TH

 [FACEBOOK.COM/PULSESCIENCEGROUP](https://www.facebook.com/pulsesciencegroup)
 [LINKEDIN.COM/COMPANY/PULSESCIENCE](https://www.linkedin.com/company/pulsescience)
 [WWW.PULSESCIENCE.CO.TH](http://www.pulsescience.co.th)

PULSE  SCIENCE

URS is a member of Registrar of Standards (Holdings) Ltd.

IncuWork *mi* is a microscope stage incubator, providing optimal stability of the environmental parameters for living cells during microscope examining

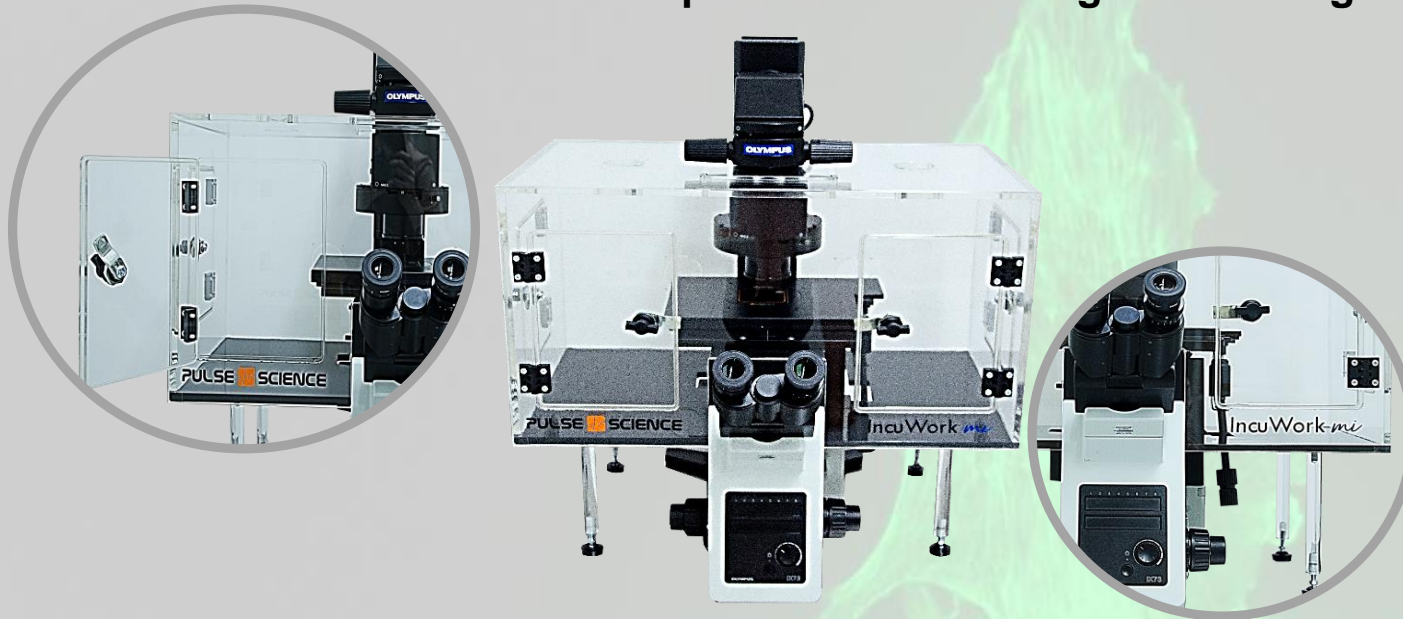

Maintaining optimal temperature for living cells on microscope stage

- Temperature Control Unit generates warm air to provide the temperature inside the chamber.
- Double warm air input ensure temperature uniformity inside the chamber
- Stable and precise temperature throughout the chamber, range from 5 above ambient to 45
- Temperature accuracy $\pm 0.1^{\circ}\text{C}$ at the sample

The ideal solution for living cell imaging

- Prepared for any inverted microscope

Ergonomic design

- The focusing knob and stage control are outside allowing easy control of microscope
- Large front and side door to ensure comfort to access sample
- The heater is placed outside ensuring minimum vibration to the microscope
- No affect on focus stability

Option feature

- Small gas chamber and humidity module can be added to ensure optimal environment for living cells.

